


# BASIC Table

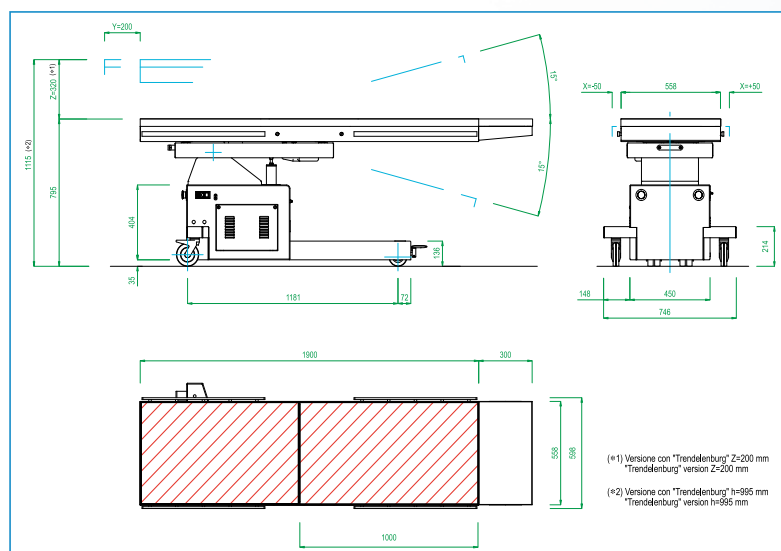

## CONTATTI

**INTERMEDICAL S.r.l.**  
IMD Group

24050 Grassobbio (Bergamo) Italy  
Via E. Fermi, 26  
ph. 0039 035 6594811  
fax 0039 035 6594899  
e-mail [info@inter-med.it](mailto:info@inter-med.it)  
[www.inter-med.it](http://www.inter-med.it)

Data may be altered with no prior warning

doc. 11425 - 27/06/2013



IMPROVING PERFORMANCES



*Basic Table is a mobile and wheeled multifunctional table. Thanks to the radio-transparent plane and the versatility of its performances, it ensures excellent results both in fluoroscopy and radiography examinations. Fit for a wide range of applications, it offers an ideal use in conjunction with C-arm Mobile Image Intensifiers and any other mobile or stationary x-ray system. Basic Table may be classified as a universal table and it is perfectly suitable to the widest range of applications, such as gynaecology, urology, lithotripsy and surgery.*

*It is provided with a remote control handswitch for the motorized movements: up/down ("Z" axis= 320 mm), trendelenburg (+/- 15°), "Y" axis movement (200 mm), "X" axis movement (+/- 50 mm). The up/down movement can be also activated by a foot-switch.*



*The cushion set is made in fire-retardant and biocompatible PVC (leather imitation), but it is available, on request, the anti-static and electrically conductive version. In order to satisfy every user's requirements, Basic Table is available in the version with up/down and X/Y movements and in the version with up/down, X/Y movements and trendelenburg.*

*Each Basic Table is also provided with steel bars for the fixing of some optional accessories, such as 300 mm table extension, leg-rest, shoulder-rest, arm-rest, phlebo stand, roll holder, head-rest and liquid tank with drain pipe. Besides, for the lithotripsy use, it is possible to install proper supports for the treatment with shock wave membrane.*

